



For Professional Use Only

# TMPRSS2 Break Apart FISH Probe Kit

#### Introduction


The TMPRSS2 Break Apart FISH Probe Kit is designed to detect rearrangements in the human *TMPRSS2* gene located on chromosome band 21q22.3. In addition to revealing breaks, which can lead to translocation of parts of the gene, inversion, or its fusion to other genes, the probe set can also be used to identify other *TMPRSS2* aberrations such as deletions or amplifications. Rearrangements and abnormal expression of the *TMPRSS2* gene – also known as *PP9284* or *PRSS10* – have been observed in prostate cancer and other malignancies.

### **Intended Use**

To detect rearrangements in the human *TMPRSS2* gene located on chromosome band 21q22.3.

Cont.	Color
LSP TMPRSS2 5' FISH Probe	CytoGreen
LSP TMPRSS2 3' FISH Probe	CytoOrange

## **Probe Design**

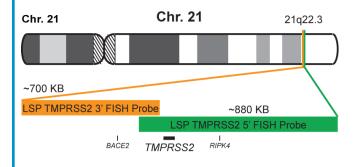

LSP TMPRSS2 5' FISH Probe covers the entire *TMPRSS2* gene, the 5' (start) portion of the gene and some adjacent genomic sequences; it also covers some sequences adjacent to the 3' portion of the gene. LSP TMPRSS2 3' FISH Probe covers sequences downstream of the 3' end of the gene. The two probes are flanking sequences across the *TMPRSS2* gene in which variable breakpoints have been observed.

Not to Scale

Cat. No.	Volume
CT-PAC154-10-GO	10 Tests (100 μL)

## **Signal Pattern Interpretation**

Normal Pattern

20G\*

Abnormal Pattern

Other Patterns

\*Overlapping orange and green signals can appear as yellow.

- 1) Paoloni-Giacobino A, et al. Genomics. 44(3):309-20 (1997).
- 2) Vaarala MH, et al. *Int J Cancer*. 94(5):705-10 (2001).
- 3) Vaarala MH, et al. *J Pathol.* 193(1):134-40 (2001).
- 4) Choi SY, et al. *Trends Mol Med.* 15(7):303-12 (2009)
- 5) Barwick BG, et al. *Br J Cancer*. 102(3):570-6 (2010).

CytoTest Inc. 9430 Key West Ave., Suite 210 Rockville, MD 20850, USA

#### Life Technologies (India) Pvt. Ltd.

306, Aggarwal City Mall, Opposite M2K Pitampura, Delhi – 110034 (INDIA). Ph: +91-11-42208000, 42208111, 42208222, Mobile: +91-9810521400, Fax: +91-11-42208444 Email: <a href="mailto:customerservice@lifetechindia.com">customerservice@lifetechindia.com</a> Website: <a href="https://www.lifetechindia.com">www.lifetechindia.com</a>

<sup>\*</sup> CE IVD only available in certain countries. All other countries are either ASR or RUO. Please contact your local dealer or our headquarters for more information.

India Contact: