



For Professional Use Only

# IGH-MYC Fusion/Translocation FISH Probe Kit

## Introduction


The IGH-MYC Fusion/Translocation FISH Probe Kit is designed to detect rearrangements involving the human IGH locus and MYC gene, located on chromosome bands 14q32.33 and 8q24.21, respectively. Rearrangements between the two regions have been observed in Burkitt's Lymphoma (BL) and other lymphomas and leukemias.

## **Intended Use**

To detect rearrangements involving the human /GH locus and /MYC gene located on chromosome bands 14q32.33 and 8q24.21, respectively.

Cont.	Color
LSP IGH 5'-3' FISH Probe	CytoGreen
LSP MYC 5'-3' FISH Probe	CytoOrange

## **Probe Design**

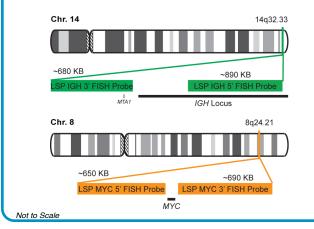

LSP IGH 5'-3' FISH Probe covers the 5' and the center sequences of the *IGH* locus, and it also covers the 3' (end) part and the neighboring downstream region. LSP MYC 5'-3' FISH Probe covers some genomic sequences downstream of the 5' (start) portion of the *MYC* gene, and it also covers the sequences adjacent to the 3' end of the gene. The probe set is optimized to reveal translocations between the two regions.

Cat. No.	Volume
CT-PAC223-10-GO	10 Tests (100 μL)

## **Signal Pattern Interpretation**

Normal Pattern
2O + 2G\*

Abnormal Pattern
Other Patterns

\*Overlapping orange and green signals can appear as yellow.

- 1) Depinho RA, et al. *Ann Clin Res.* 18(5-6):284-9 (1986).
- 2) Garte SJ. Crit Rev Oncog. 4(4):435-49 (1993).
- 3) Einerson RR, et al. *Leukemia*. 20:1790–9 (2006)
- 4) Le Gouill S, et al. *Haematologica*. 92(10):1335-42 (2007).
- 5) Blancato J, et al. *Br J Cancer*. 90(8):1612-9 (2004).

CytoTest Inc. 9430 Key West Ave., Suite 210 Rockville, MD 20850, USA

## India Contact:

#### Life Technologies (India) Pvt. Ltd.

<sup>\*</sup> CE IVD only available in certain countries. All other countries are either ASR or RUO. Please contact your local dealer or our headquarters for more information.