



Protein G Magnetic Beads comparison

1. Abraxis AbraMag Protein G Beads
2. Competitor 1: Bangs Laboratories Protein G BioMag particles
3. Competitor 2: ThermoFisher Pierce Protein G beads

Antibody: Rabbit IgG 10.7mg/mL, Sigma; 12.5µg, 25µg, 50µg added
0.1 mg beads used for all samples

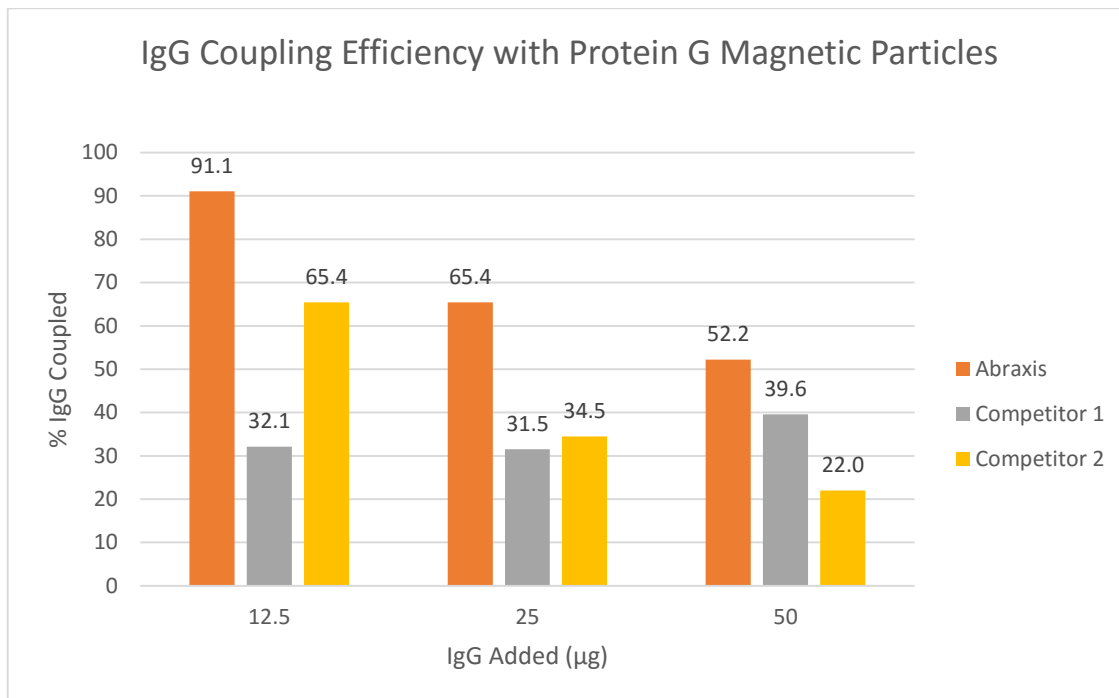

SDS-PAGE of Flow-through

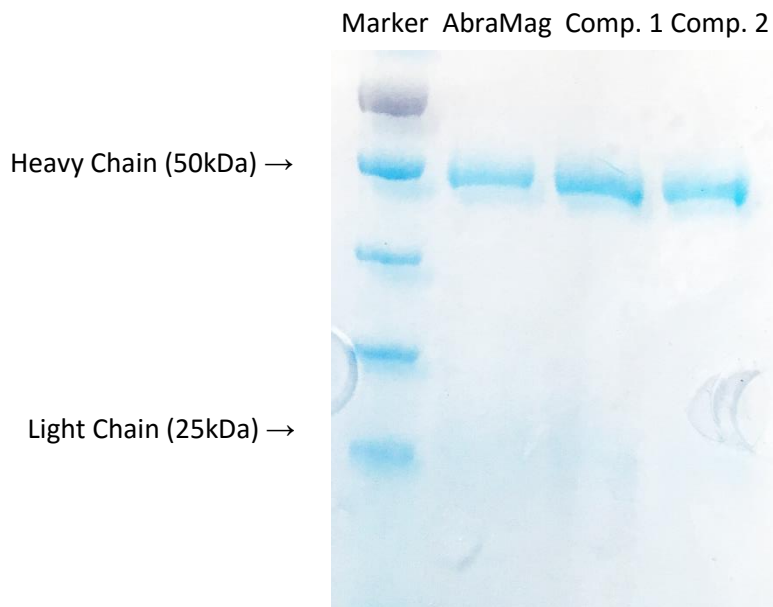

Amount of non-bound antibody after adding 0.5 mg of rabbit IgG to 1 mg Protein G beads.

Note that less antibody is lost using AbraMag beads compared to competitor beads.

India Contact:

Life Technologies (India) Pvt. Ltd.

306, Aggarwal City Mall, Road No. 44, Pitampura, Delhi – 110034, India

Mobile: +91-98105-21400, Tel: +91-11-42208000, 8111, 8222, Fax: +91-11-42208444

Email: customerservice@lifetechindia.com, www.atzlabs.com ; www.lifetechindia.com